

Prisms & Polarizers

Prisms

| | |
|--|----|
| <i>Introduction to Right Angle Prisms</i> | 72 |
| <i>Laser Quality Right Angle Prisms</i> | 73 |
| <i>Laboratory Quality Right Angle Prisms</i> | 73 |
| <i>Equilateral Prisms</i> | 74 |
| <i>Wedge Prisms</i> | 75 |
| <i>Anamorphic Prism Pair</i> | 75 |
| <i>Penta Prisms</i> | 76 |
| <i>Dove Prisms</i> | 76 |
| <i>Brewster Prisms</i> | 77 |
| <i>Littrow Prisms</i> | 77 |

Polarizers

| | | |
|---------------------------|--|----|
| Laser Polarizers | <i>Broadband Polarizing Beamsplitting Cubes</i> | 78 |
| | <i>Narrowband Polarizing Beamsplitting Cubes</i> | 79 |
| Calcite Polarizers | <i>Introduction to Calcite Polarizers</i> | 80 |
| | <i>Glan Laser Calcite Polarizers</i> | 80 |
| | <i>High Transmission Glan Laser Calcite Polarizers</i> | 81 |
| | <i>Glan Taylor Calcite Polarizers</i> | 81 |
| | <i>Glan Thompson Calcite Polarizers</i> | 82 |
| | <i>Glan Thompson Calcite Beamsplitting Prisms</i> | 82 |
| | <i>Soleil-Babinet Compensator</i> | 83 |
| Sheet Polarizers | <i>Introduction to Sheet Polarizers</i> | 84 |
| | <i>UV-Visible Linear Polarizers</i> | 84 |
| | <i>Visible Linear Polarizers</i> | 85 |
| Waveplates | <i>Introduction to Waveplates</i> | 86 |
| | <i>Zero Order Waveplates</i> | 87 |
| | <i>Multiple order Waveplates</i> | 87 |

Optics

Lenses & Microscope Components

Coatings

Mirrors, Beamsplitters & Windows

Prisms & Polarizers

Filters

Pinholes

Opto-mechanics

Tables, Breadboards & Rails

Mounting Hardware

Mirror & Component Mounts

Manual Micro-positioners

Motorized Positioners

Lasers & Accessories

Beam Delivery


Laser Measurement

Diode Laser Modules

Broadband Polarizing Beamsplitter Cubes



- Laser quality
- Usable over wide wavelength ranges

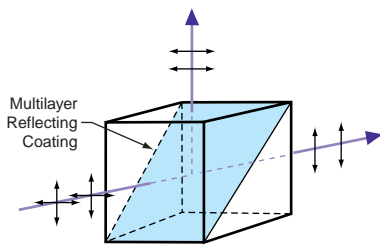

Specifications

Material: SF-2 glass
Transmission (p-polarized): >90% average
Reflection (s-polarized): >99.8% average
AR coating: R < 0.5% per surface
Transmitted wavefront: $\lambda/4$ at 633 nm
Surface quality: 20-10
Extinction ratio: >500:1
Dimensions: ± 0.2 mm
Clear aperture: 85% of cube dimension
Laser damage threshold
CW: 2 kW/cm²
Pulsed (10 ns): 1 J/cm²


These cemented Polarizing Beamsplitter Cubes are coated to enable operation over a wide range of wavelengths.

The polarization separation is excellent with the transmitted and reflected beams at 90° to each other irrespective of wavelength.

Broadband Polarizing Beamsplitter Cubes

| Wavelength Range (nm) | 12.7 mm Cube | | 25.4 mm Cube | |
|-----------------------|----------------|----------|----------------|----------|
| | Catalog Number | Price US | Catalog Number | Price US |
| 450-700 | 44-4703 | \$231.00 | 44-4711 | \$310.00 |
| 670-980 | 44-4729 | \$242.00 | 44-4737 | \$320.00 |
| 1300-1550 | 44-4745 | \$242.00 | 44-4752 | \$320.00 |

Broadband Laser Optics

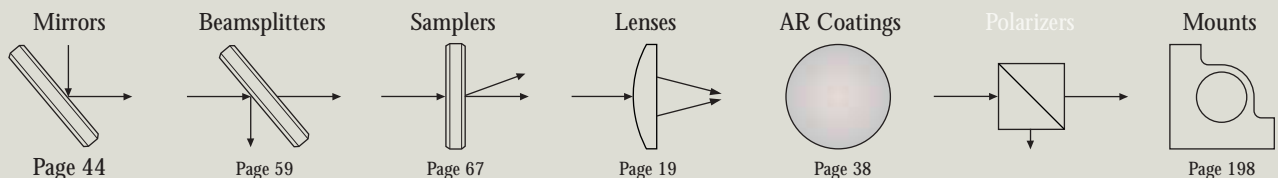

Optics

Lenses & Microscope Components

Coatings

Mirrors, Beamsplitters & Windows

Prisms & Polarizers

Filters

Pinholes

Opto-mechanics

Tables, Breadboards & Rails

Mounting Hardware

Mirror & Component Mounts

Manual Micro-positioners

Motorized Positioners

Lasers & Accessories


Beam Delivery

Laser Measurement

Diode Laser Modules


Narrowband Polarizing Beamsplitter Cubes

- Laser quality
- Common laser wavelengths

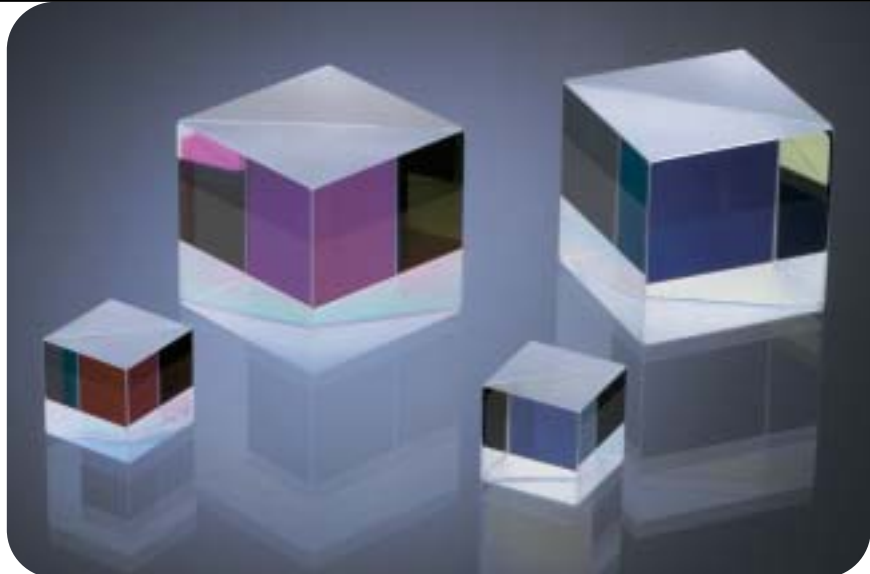

- Optics
- Lenses & Microscope Components
- Coatings
- Mirrors, Beamsplitters & Windows
- Prisms & Polarizers
- Filters
- Pinholes

Narrowband polarizing beamsplitter cubes split the polarization states of an incoming beam. The P-polarization passes straight through the cube, while the S-polarization is reflected 90°. This split is the result of the multi-layer dielectric coating that has been placed on the internal hypotenuse surface. The entrance and exit faces have been AR coated with multi-layer coatings to maximize transmission.

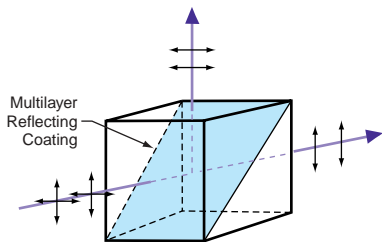

Kinematic Prism Platforms



Kinematic Prism Platforms are also available.

- Opto-mechanics
- Tables, Breadboards & Rails
- Mounting Hardware
- Mirror & Component Mounts
- Manual Micro-positioners
- Motorized Positioners

Specifications

Material: BK7 glass
Transmission (p-polarized): >95%
Reflection (s-polarized): >99.8%
AR Coating: R ≤ 0.25% per surface
Typical Useable Bandwidth: 10 nm
Transmitted Wavefront: λ/4 at 633nm
Surface Quality: 20-10
Extinction Ratio: 1000:1
Dimensions: ±0.2 mm
Clear Aperture: 80% of cube dimension
Laser Damage Threshold
CW: 2 kW/cm²
Pulsed (10 ns): 1 J/cm²

Narrowband Polarizing Beamsplitter Cubes

| Wavelength (nm) | 12.7 mm Cube | | 25.4 mm Cube | |
|-----------------|----------------|----------|----------------|----------|
| | Catalog Number | Price US | Catalog Number | Price US |
| 488 | 44-4380 | \$190.00 | 44-4398 | \$270.00 |
| 514 | 44-4406 | \$190.00 | 44-4414 | \$270.00 |
| 532 | 44-4422 | \$190.00 | 44-4430 | \$270.00 |
| 633 | 44-4448 | \$190.00 | 44-4455 | \$270.00 |
| 670 | 44-4489 | \$190.00 | 44-4497 | \$270.00 |
| 780 | 44-4505 | \$190.00 | 44-4513 | \$270.00 |
| 808 | 44-4521 | \$190.00 | 44-4539 | \$270.00 |
| 830 | 44-4547 | \$190.00 | 44-4554 | \$270.00 |
| 850 | 44-4562 | \$190.00 | 44-4570 | \$270.00 |
| 1064 | 44-4588 | \$190.00 | 44-4596 | \$270.00 |
| 1300 | 44-4604 | \$190.00 | 44-4612 | \$270.00 |
| 1550 | 44-4620 | \$190.00 | 44-4638 | \$270.00 |

- Lasers & Accessories
- Beam Delivery
- Laser Measurement
- Diode Laser Modules

Calcite Polarizers

Optics

Lenses & Microscope Components

Coatings

Mirrors, Beamsplitters & Windows

Prisms & Polarizers

Filters

Pinholes

- **Broadband linear polarizers**
- **High extinction ratio**
- **Laser quality**

Calcite Polarizers are laser quality, linear polarizers. They can be used over a broad wavelength region, making them suitable for multiple applications. Calcite Polarizers consist of two prisms made of a high birefringent calcite material. The internal interface between the prisms is cut at Brewster's angle. This design allows for linear polarization of the incoming light with a high extinction ratio over a very wide wavelength region.

Ealing offers five different types of calcite polarizers, each with their own unique advantages. Care should be taken to select the best one for the application at hand based on the lasers' power density, wavelength, and required acceptance angle. All Calcite polarizers are supplied cemented in a slotted black anodized aluminum cylinder, which can be mounted with polarizers mounts.

Opto-mechanics

Tables, Breadboards & Rails

Mounting Hardware


Mirror & Component Mounts

Manual Micro-positioners

Motorized Positioners

Glan Laser Polarizers

- **Ideal for high power lasers**
- **Low beam divergence**

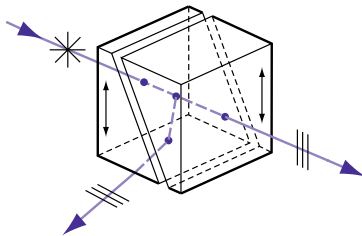

Glan Laser polarizers are intended for use with high power laser beams. Their many applications include intra-cavity gain switching, beam combination, pulse extraction and feedback elimination. Glan Laser Polarizers are air spaced, and have side windows to enable rejected beams to escape, thus avoiding damage due to excess heating. Care must be taken to terminate the rejected beams, though.

Glan Laser Polarizers


| Catalog Number | Clear Aperture (mm) | Price US |
|----------------|---------------------|------------|
| 43-6741 | 10.0 | \$775.00 |
| 43-6766 | 15.0 | \$1,175.00 |
| 43-6774 | 20.0 | \$1,935.00 |

Specifications

Material: Laser grade calcite
Wavelength Range: 220-2800 nm
Peak Transmission: 88%
Extinction Ratio: 10^5
Beam Deviation: <3 arcmin
Flatness: $\lambda/8$ at 589 nm
Surface Quality: 20-10
Field Angle: 8°

| Mount Dimensions | Diameter (mm) | Length (mm) |
|------------------------|---------------|-------------|
| 10 mm Aperture: | 25.4 | 32.0 |
| 15 mm Aperture: | 31.8 | 38.0 |
| 20 mm Aperture: | 38.1 | 44.0 |

Dimensions Tolerance: ± 0.1 mm
Laser Damage Threshold
CW: 100 W/cm²
Pulsed (1ns): 300 MW/cm²

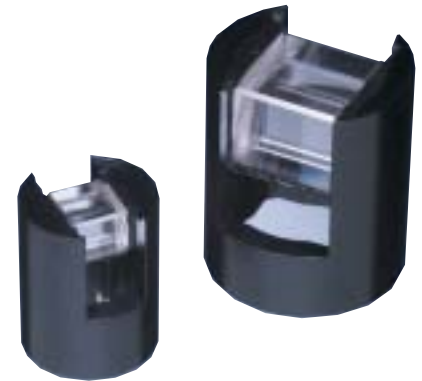

Calcite Polarizer Mounts



Calcite Polarizer Mounts are also available.

Lasers & Accessories


Beam Delivery

Laser Measurement

Diode Laser Modules

High Transmission Glan Laser Polarizers

- Ideal for high power lasers
- High transmission

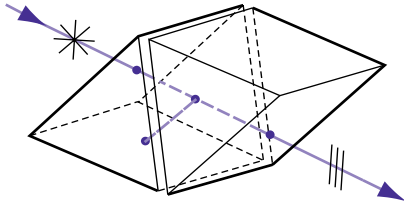

These Glan Laser polarizers are designed as high power, high transmission linear polarizers. All entrance and exit surfaces are at Brewster's angle, resulting in a maximum transmission of the p-polarized incoming light. Note that although the High Transmission Glan Laser Polarizer will work over an 8° field angle, in order to achieve 98% transmission it is necessary to strike the entrance surface at Brewster's angle.

Specifications

Material: Laser grade calcite
Wavelength range: 220-2800 nm
P-polarized Transmission: 98%
Extinction Ratio: 10⁵
Flatness: λ/8 at 589 nm
Surface Quality: 20-10
Field Angle: 8°

Beam Displacement

10mm Aperture: 5 mm
15mm Aperture: 7.5 mm

| Mount Dimensions | Diameter (mm) | Length (mm) |
|------------------|---------------|-------------|
| 10mm Aperture: | 25.4 | 32.0 |
| 12mm Aperture: | 25.4 | 34.0 |
| 15mm Aperture: | 31.8 | 38.0 |

Dimensions Tolerance: ±0.1mm

Laser Damage Threshold

CW: 100 W/cm²
Pulsed (1ns): 300 MW/cm²




High Transmission Glan Laser Polarizers

| Catalog Number | Clear Aperture (mm) | Price US |
|----------------|---------------------|------------|
| 43-6782 | 10.0 | \$936.00 |
| 43-6808 | 15.0 | \$1,365.00 |

Glan Taylor Polarizers

- Useful with low or medium power lasers

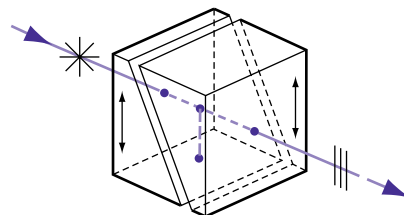

These Glan Taylor Polarizers are ideal for low and medium power applications where the rejected beam is not required. They do not have any escape windows and are assembled with fine ground, black glass cemented to the calcite prisms for efficient absorption of the rejected beam.

Specifications

Material: Optical grade calcite
Wavelength range: 220-2800 nm
Peak transmission: 88%
Extinction ratio: 10⁵
Beam deviation: <3 arcmin
Flatness: λ/8 at 589 nm
Surface quality: 20-10
Field angle: 8°

| Mount Dimensions | Diameter (mm) | Length (mm) |
|------------------|---------------|-------------|
| 10mm Aperture: | 25.4 | 32.0 |
| 15mm Aperture: | 31.8 | 38.0 |
| 20mm Aperture: | 38.1 | 44.0 |

Dimensions Tolerance: ±0.1mm

Laser Damage Threshold

CW: 10 W/cm²
Pulsed (1ns): 20 MW/cm²



Glan Taylor Polarizers

| Catalog Number | Clear Aperture (mm) | Price US |
|----------------|---------------------|------------|
| 43-6824 | 10.0 | \$695.00 |
| 43-6840 | 15.0 | \$1,095.00 |
| 43-6857 | 20.0 | \$1,859.00 |

- Optics
- Lenses & Microscope Components
- Coatings
- Mirrors, Beamsplitters & Windows
- Prisms & Polarizers
- Filters
- Pinholes
- Opto-mechanics
- Tables, Breadboards & Rails
- Mounting Hardware
- Mirror & Component Mounts
- Manual Micro-positioners
- Motorized Positioners
- Lasers & Accessories
- Beam Delivery
- Laser Measurement
- Diode Laser Modules

Optics

Lenses & Microscope Components

Coatings

Mirrors, Beamsplitters & Windows


Prisms & Polarizers

Filters

Pinholes

Glan Thompson Polarizers

- Wide acceptance angles
- Low power applications

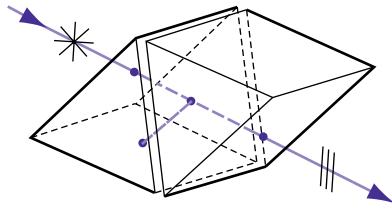

Glan Thompson Polarizers can be used over wide acceptance angles. They are available with field angles of either 15° or 26°. This makes them ideal for use with diverging or converging beams. They can be used with lower power lasers or broadband light sources.


The two calcite prisms are cemented together with index matching cement. The side faces are unpolished, and covered with a black paint to absorb the reflected component.

Specifications

Material: Optical grade calcite
Wavelength Range: 350-2500 nm
Peak Transmission: 90%
Extinction Ratio: 10⁵
Beam Deviation: <3 arcmin
Flatness: λ/8 at 589 nm
Surface Quality: 20-10

| Catalog Number | Diameter (mm) | Length (mm) |
|----------------|---------------|-------------|
| 43-6881 | 25.4 | 27.0 |
| 43-6899 | 25.4 | 32.0 |
| 43-6923 | 31.8 | 40.0 |
| 43-6931 | 31.8 | 52.0 |
| 43-6949 | 38.1 | 52.0 |

Dimensions Tolerance: ±0.1mm
Laser Damage Threshold: 1 W/cm²

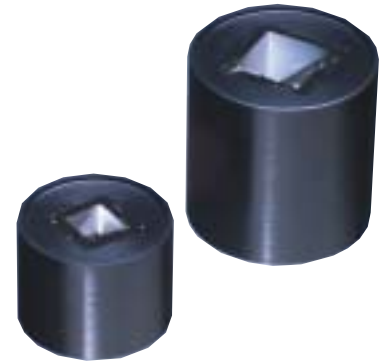

Glan Thompson Polarizers

| Catalog Number | Clear Aperture (mm) | Field Angle (deg.) | Price US |
|----------------|---------------------|--------------------|------------|
| 43-6881 | 10.0 | 15 | \$535.00 |
| 43-6899 | 10.0 | 26 | \$635.00 |
| 43-6923 | 15.0 | 15 | \$1,350.00 |
| 43-6931 | 15.0 | 26 | \$1,895.00 |
| 43-6949 | 20.0 | 15 | \$2,321.00 |

Opto-mechanics

Tables, Breadboards & Rails

Mounting Hardware

Mirror & Component Mounts


Manual Micro-positioners

Motorized Positioners

Glan Thompson Beamsplitting Prisms

Glan Thompson Beamsplitting Prisms are ideal in applications where it is essential that the transmitted p-polarized light is undeviated in its path. Unlike standard Glan Thompson Polarizers where the s-polarized ordinary ray is reflected and absorbed, these Beamsplitting Prisms have an additional

escape window to allow transmission of the ordinary ray. The escape window is designed such that the beam emerges normal to it ensuring that there is no chromatic dispersion. The p-polarized extraordinary ray is transmitted undeviated from its original path. Ealing Glan Thompson Beamsplitting Prisms have an angular deviation between the two beams of 44° which is not wavelength dependent.

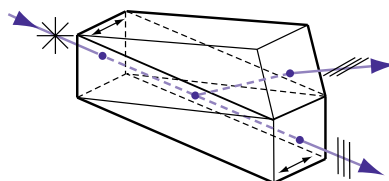

Glan Thompson Beamsplitting Prisms

| Catalog Number | Aperture (mm) | Diameter (mm) | Length (mm) | Price US |
|----------------|---------------|---------------|-------------|------------|
| 43-8523 | 10.0 | 31.8 | 48 | \$855.00 |
| 43-8531 | 12.0 | 38.1 | 49 | \$1,070.00 |

Specifications

Material: Optical calcite
Wavelength Range: 350-2500 nm
Peak Transmission: 90%
Extinction Ratio: 10⁵
Surface Quality: 20-10
Beam Deviation: <3 mins
Dimensions: ±0.1mm
Laser Damage Threshold: 1 W/cm²

Lasers & Accessories

Beam Delivery

Laser Measurement

Diode Laser Modules

Soleil-Babinet Compensator, Adjustable Retarder

- An instrument for precise determination of retardation



Specifications

Materials

Optics: Laser quality crystalline quartz
Housing: Non-magnetic metal

Diameter: 152.4 mm

Clear Aperture: 10 mm

Thickness: 105.3 mm
 (compensator and mount)

Dimensions Tolerance: ± 0.25 mm

Optical Center: 75.2 mm
 (bottom to optical center line)

Maximum Height: 243.5 mm
 (micrometer fully extended)

Wavefront Distortion: $< \lambda/4$ at 633 nm

Wavelength Range: 200-2700 nm

Retardation Range:

320 nm = 0 to 2λ
 633 nm = 0 to 1λ
 1230 nm = 0 to $\lambda/2$
 2500 nm = 0 to $\lambda/4$

Readout: 4 place digital;
 fifth place interpolated

Resolution: 0.001 λ
 (0 to 25,000 on digital readout)

Temperature Range: 0°C to 70°C

A Soleil-Babinet Compensator functions as an adjustable zero order retarder over the wavelength range 200-2700 nm. It allows complete analysis and selection of the state of polarization of a beam and can be used for inspection and comparative work.

The Ealing Soleil-Babinet Compensator consists of two crystalline quartz wedges with their optic axes parallel and at 45° to the polarization direction of the input beam. One wedge is fixed and the other, which is attached to a crystalline quartz compensating block with its axis at 90° to the wedge, is adjustable by a micrometer screw. This adjustment changes the path difference through the instrument and hence the retardation.

It can be used to select a uniform phase difference between the extraordinary and ordinary rays of 0- 2π . As a result incident elliptically or circularly polarized light can be converted into linearly polarized light by introducing the appropriate

compensation. Conversely any desired polarization form can be obtained by pre-setting the appropriate values.

The Soleil-Babinet Compensator is mounted conveniently on a precision ballbearing indexing head which has a fixed outer circumference graduated 0°, 180°, +45°, +90°, +135°, -45°, -90° and -135°. The inner circumference is rotatable through 360° and has indicator marks at one degree increments with each 10° being labelled. The outer circumference has a knurled locking screw for absolute fixing. A 1/4-20 tapped hole is located at the 180° mark for post mounting. The micrometer adjustment screw has a four-place digital readout. The fifth place can be interpolated from alignment marks on the micrometer barrel. All Compensators are supplied in wooden instrument cases with instruction manual and calibration data.

Soleil-Babinet Compensator — Adjustable Retarder

| Catalog Number | Price US |
|----------------|------------|
| 34-5918 | \$3,450.00 |

Optics

Lenses & Microscope Components

Coatings

Mirrors, Beamsplitters & Windows

Prisms & Polarizers

Filters

Pinholes

Opto-mechanics

Tables, Breadboards & Rails

Mounting Hardware

Mirror & Component Mounts

Manual Micro-positioners

Motorized Positioners

Lasers & Accessories

Beam Delivery

Laser Measurement

Diode Laser Modules

Sheet Polarizers

Optics

Lenses & Microscope Components

Coatings

Mirrors, Beamsplitters & Windows

Prisms & Polarizers

Filters

Pinholes

Opto-mechanics

Tables, Breadboards & Rails

Mounting Hardware

Mirror & Component Mounts

Manual Micro-positioners

Motorized Positioners

Lasers & Accessories

Beam Delivery

Laser Measurement

Diode Laser Modules

- *Ideal for broadband applications*
- *Cost effective*

For low power applications Sheet Polarizers can often provide a simple cost effective solution.

Unpolarized light passing through a Sheet Polarizer emerges as linearly polarized light. Sheet Polarizers are often used in pairs. When the two polarizers have their axes aligned transmission is at a maximum. When the axes are orthogonal to each other transmission is nearly zero. For intermediate positions the transmission is given by the equation:


$$I = I_{\max} \cos^2\theta$$

where θ = angle between the axes of the polarizers.

UV-Visible Linear Polarizers

- *Useful for 200-800 nm*
- *UV fused silica substrate*

Providing excellent linear polarization from 200-800 nm in a single special coating, these polarizers are mounted on UV quality, fused silica. The coated surface is delicate and care should be taken while handling or cleaning. Normal cleaning solvents should not be used. For best results only clean with blown air. The coating is neutral green in color. Maximum survival temperature is 95°C. Polarization and transmission characteristics do not

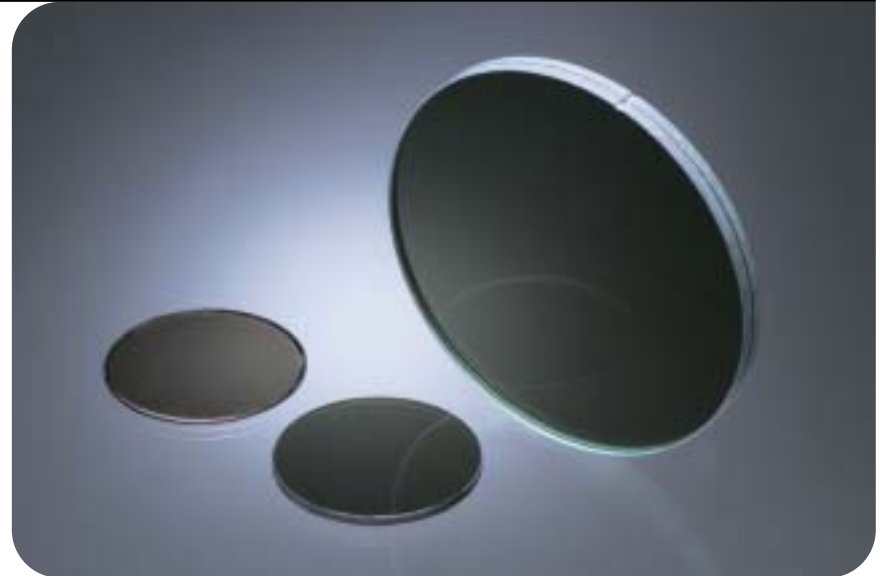

Ealing offers Sheet Polarizers for the visible, ultraviolet and infrared spectrums.


Sheet Polarizers are not suitable for high power use. Their main applications are for production and analysis of polarized radiation, and the elimination of unwanted beams.

The degree of polarization is virtually independent of the incident angle. This allows them to be used with highly convergent or divergent beams and still produce uniform polarization.

deteriorate with time or long exposure to UV, visible, or IR irradiation. Thickness is nominally 1.6 mm.

UV-Visible Linear Polarizers

| Catalog Number | Diameter (mm) | Price US |
|----------------|---------------|------------|
| 23-2520 | 25.4 | \$407.00 |
| 23-2363 | 50.8 | \$620.00 |
| 22-9039 | 101.6 | \$1,554.00 |

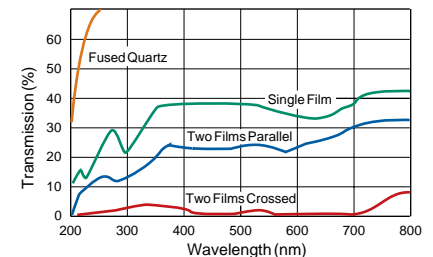


Visible Linear Polarizers

- Usable over 350-750 nm
- Ideal for broadband applications


These sheet polarizers are designed for 350-750 nm. Total luminous transmittance is 22% for white light; total integrated transmission for two crossed polarizers is 0.05%. The polarizing sheet is sandwiched between two optically ground and polished pieces of glass. Polarization and transmission characteristics are stable with time and prolonged visible

irradiation. High intensity ultraviolet or infrared irradiation will degrade both the polarizing properties and total transmission of the unit. Recommended temperature extremes are -60°C and +80°C. Maximum survival temperature is 90°C. High relative humidity will tend to cause a separation of the glass sandwich and should be avoided.

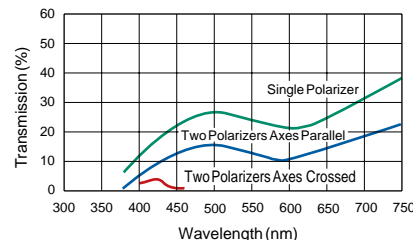

Visible Linear Polarizers

| Catalog Number | Diameter (mm) | Thickness (mm) | Price US |
|----------------|---------------|----------------|----------|
| 23-5671 | 50.8 | 3.2 | \$280.00 |
| 22-9062 | 101.6 | 7.2 | \$357.00 |

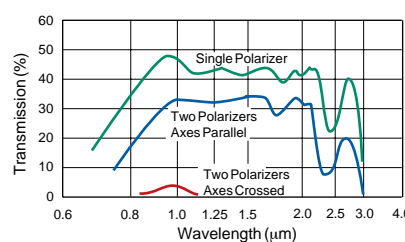

Rotatable Mounts



Optics

Lenses & Microscope Components

Coatings

Mirrors, Beamsplitters & Windows

Prisms & Polarizers

Filters

Pinholes

Opto-mechanics

Tables, Breadboards & Rails

Mounting Hardware

Mirror & Component Mounts

Manual Micro-positioners

Motorized Positioners

Lasers & Accessories

Beam Delivery

Laser Measurement

Diode Laser Modules

Waveplates

Optics

Lenses & Microscope Components

Coatings

Mirrors, Beamsplitters & Windows

Prisms & Polarizers

Filters

Pinholes

Opto-mechanics

Tables, Breadboards & Rails

Mounting Hardware

Mirror & Component Mounts

Manual Micro-positioners

Motorized Positioners

Lasers & Accessories

Beam Delivery

Laser Measurement

Diode Laser Modules

- Quarter and half waveplates
- Zero or multiple order available

Applications for Half Waveplates include rotating the plane of polarization (e.g. in a laser), electro-optic modulation and as a variable ratio beamsplitter (when used in conjunction with a polarizing cube).

Waveplates are made from materials which exhibit birefringence. The velocities of the extraordinary and ordinary rays through the birefringent material vary inversely with their refractive indices. For the case of crystal quartz the extraordinary beam has a higher refractive index and therefore a slower velocity. For this reason its direction is known as the 'slow' axis. Likewise, the direction of the ordinary beam is known as the 'fast' axis and is indicated by the marked lines on the mount.

The difference in velocities gives rise to a phase difference when the two beams recombine. In the case of an incident linearly polarized beam this is given by

$$\theta = \frac{\pm 2\pi d (n_e - n_o)}{\lambda}$$

where:

θ = phase difference

d = thickness of waveplate in mm


n_e, n_o = refractive indices of extraordinary and ordinary rays respectively

λ = wavelength in nm

At any specific wavelength the phase difference is governed by the thickness of the retarder. Quarter and Half Waveplates are two specific cases of this.



Quarter Waveplate

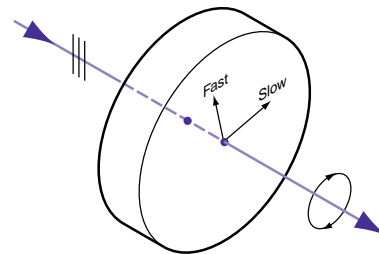


A Quarter Waveplate is used to convert linearly polarized beams into circularly polarized beams (and vice versa).

The construction of a Quarter Waveplate is such that the fast axis, indicated by the marked line, lies in the surface at 45° to the input polarization. The input beam is resolved into two components of equal amplitude, but each with a different velocity.

Applications for Quarter Waveplates include creating circular polarization from linear or linear polarization from circular, ellipsometry, optical pumping, suppressing unwanted reflections (when used in conjunction with a polarizer) and optical isolation (when used with a Polarizing Beamsplitter Cube.)

Half Waveplates

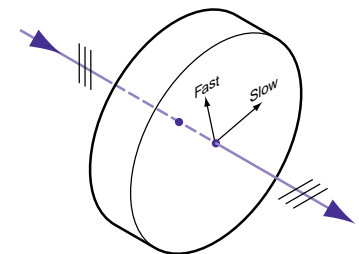

The thickness of a Half Waveplate is such that the phase difference is π (zero order) or 3π , 5π , 7π , etc (multiple orders). A linearly polarized beam incident on a Half Waveplate emerges as a linearly polarized beam but rotated such that its angle to the optic axis is twice that of the incident beam. It is usual to have the fast axis lying in the surface of the retarder at 45° to the input polarization. The Half Waveplate therefore introduces a 90° rotation of the plane of polarization.

Zero Order Waveplates

Ealing Zero Order Waveplates are the preferred type of waveplate. They are not as sensitive to changes in temperature, wavelength, angle of incidence or collimation. A wavelength shift of 15 nm will result in approximately a 1% retardation change. They are supplied in a 25.4 mm mount.

Specifications

Material: Crystalline quartz

Retardation: $\pm 0.005\lambda$

Wavefront Distortion: $\lambda/8$

AR Coating: $<0.25\%R$ per surface

Diameter: 25.4 +0.0 -0.15 mm

Aperture: 15 mm

Thickness: 8 +0.0 -0.25 mm

Laser Damage Threshold: 500 MW/cm²

Zero Order Waveplates

| Wavelength (nm) | Quarter Waveplate | | Half Waveplate | |
|-----------------|-------------------|----------|----------------|----------|
| | Catalog Number | Price US | Catalog Number | Price US |
| 248 | 45-7564 | \$372.00 | 45-7788 | \$372.00 |
| 266 | 45-7572 | \$372.00 | 45-7796 | \$372.00 |
| 308 | 45-7580 | \$372.00 | 45-7804 | \$372.00 |
| 355 | 45-7598 | \$372.00 | 45-7812 | \$372.00 |
| 488 | 45-7606 | \$372.00 | 45-7820 | \$372.00 |
| 514 | 45-7614 | \$372.00 | 45-7838 | \$372.00 |
| 532 | 45-7622 | \$372.00 | 45-7846 | \$372.00 |
| 633 | 45-7648 | \$372.00 | 45-7861 | \$372.00 |
| 670 | 45-7663 | \$372.00 | 45-7887 | \$372.00 |
| 780 | 45-7697 | \$372.00 | 45-7911 | \$372.00 |
| 808 | 45-7705 | \$372.00 | 45-7929 | \$372.00 |
| 820 | 45-7713 | \$372.00 | 45-7937 | \$372.00 |
| 830 | 45-7721 | \$372.00 | 45-7945 | \$372.00 |
| 850 | 45-7739 | \$372.00 | 45-7952 | \$372.00 |
| 905 | 45-7747 | \$372.00 | 45-7960 | \$372.00 |
| 1064 | 45-7754 | \$372.00 | 45-7978 | \$372.00 |
| 1300 | 45-7762 | \$372.00 | 45-7986 | \$372.00 |
| 1550 | 45-7770 | \$372.00 | 45-7994 | \$372.00 |

Multiple Order Waveplates

Ealing Multiple Order Waveplates are available in a range of laser-line wavelengths, mounted in a 25.4 mm diameter mount. Multiple Order Waveplates are sensitive to changes in temperature, angle or incidence and degree of collimation. They are intended for use at the design wavelength only. If used at a different wavelength, a retardation change of 10% occurs for every 0.2 nm deviation.

Specifications

Material: Crystalline quartz

Retardation: $\pm 0.005\lambda$

Wavefront Distortion: $\lambda/8$

AR Coating: $<0.25\%R$ per surface

Diameter: 25.4 +0.0 -0.15 mm

Aperture: 15 mm

Thickness: 8 +0.0 -0.25 mm

Laser Damage Threshold: 500 MW/cm²

Multiple Order Waveplates

| Wavelength (nm) | Quarter Waveplate | | Half Waveplate | |
|-----------------|-------------------|----------|----------------|----------|
| | Catalog Number | Price US | Catalog Number | Price US |
| 248 | 45-6848 | \$305.00 | 45-7200 | \$305.00 |
| 266 | 45-6855 | \$305.00 | 45-7069 | \$305.00 |
| 308 | 45-6863 | \$305.00 | 45-7077 | \$305.00 |
| 355 | 45-6871 | \$305.00 | 45-7085 | \$305.00 |
| 488 | 45-6889 | \$305.00 | 45-7093 | \$305.00 |
| 514 | 45-6897 | \$305.00 | 45-7101 | \$305.00 |
| 532 | 45-6905 | \$305.00 | 45-7119 | \$305.00 |
| 633 | 45-6921 | \$305.00 | 45-7135 | \$305.00 |
| 670 | 45-6947 | \$305.00 | 45-7150 | \$305.00 |
| 780 | 45-6970 | \$305.00 | 45-7184 | \$305.00 |
| 808 | 45-6988 | \$305.00 | 45-7192 | \$305.00 |
| 820 | 45-7002 | \$305.00 | 45-7218 | \$305.00 |
| 830 | 45-7010 | \$305.00 | 45-7226 | \$305.00 |
| 850 | 45-7016 | \$305.00 | 45-7234 | \$305.00 |
| 905 | 45-7044 | \$305.00 | 45-7259 | \$305.00 |
| 1064 | 45-7028 | \$305.00 | 45-7238 | \$305.00 |
| 1300 | 45-7036 | \$305.00 | 45-7242 | \$305.00 |
| 1550 | 45-7048 | \$305.00 | 45-7262 | \$305.00 |

Optics

Lenses & Microscope Components

Coatings

Mirrors, Beamsplitters & Windows

Prisms & Polarizers

Filters

Pinholes

Opto-mechanics

Tables, Breadboards & Rails

Mounting Hardware

Mirror & Component Mounts

Manual Micro-positioners

Motorized Positioners

Lasers & Accessories

Beam Delivery

Laser Measurement

Diode Laser Modules